


# NOW MAKING AND MOUNTING BEAM SHAPING BLOCKS IS EVEN EASIER

## *APP-1* AUTOMATIC SHADOW TRAY PLOTTER

Draw shielding block outlines directly on shadow tray with exclusive, specially designed, plotter and software.

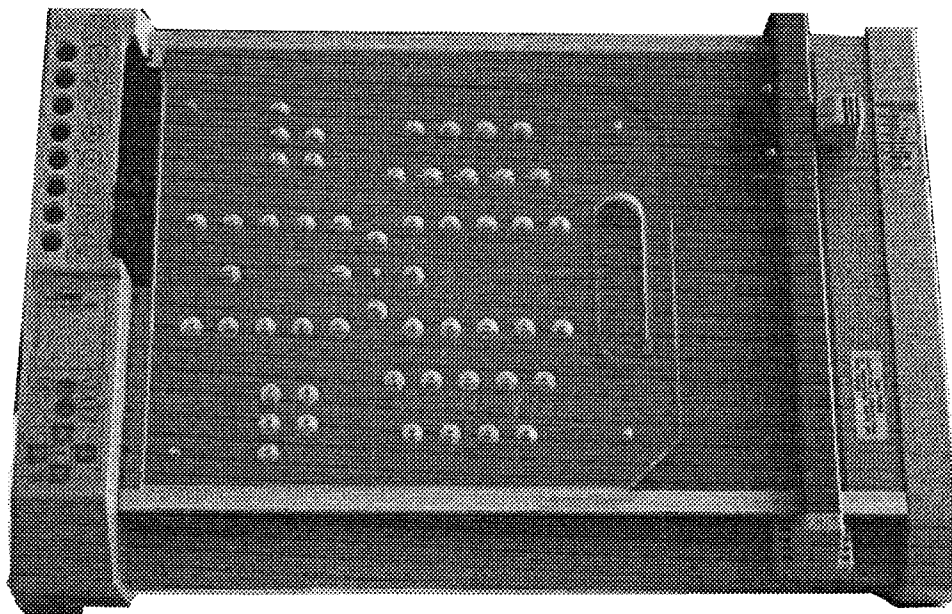

*Accurate, fast positioning of shielding blocks on a tray plate.*

The Par Scientific ACD-4MK3 and Dicut systems will automatically cut the block shapes from the digitized or planning system data and the NEW APP-1 plate plotter will draw the block outlines and the boundary of the field size directly on the plastic tray plate for fast precise mounting.

Par Scientific custom builds the plotter to fit the thickness of the tray plates used by your department!

The DOS and Windows versions of our software enable you to have up to 10 different user defined center settings for 10 different tray plate sizes. By digitizing only two extra points, you can have the actual treatment field size enclosed. The correctly scaled field size frame can be plotted on a tray plate or on paper.

The **ACD-4**, **Dicut** and **APP-1 Plotter** utilize the latest developments in digitization and automation and eliminate one of the last remaining manual tasks in the modern high technology treatment center.

S&S Par Scientific Inc.  
1101 Linwood Street  
Brooklyn N.Y. 11208  
[www.ssxray.com](http://www.ssxray.com)  
800/347-XRAY  
718/649-8500 Fax 718/257-0219

